

āfex Active Starter Kit

Available on prescription from your GP or by ordering privately on www.iMEDicare.co.uk or calling +441923237795

afex Active – Contents:	Variant – Brief Size (above Hip)	NHS Code:
1 x High Receptacle	Small (28-30")	A100-S-08
1 x Briefs (Variable Size)	Medium (32-34")	A100-M-08
1 x Re-useable 500 ml Collection Bag	Large (36-38")	A100-L-08
1 x Cleansing Container	XL (40-42")	A100-XL-08
1 x Cleansing Solution – 118ml	2X (44-46")	A100-2X-08
1 x Instructional DVD	3X (48-50")	A100-3X-08



Black Brief option prescribable separately

afex Core Supporter

Starter Kit (high pubic pressure)

Available on prescription from your GP or by ordering privately on www.iMEDicare.co.uk or calling +441923237795

afex Core Supporter – Contents:	Variant – Brief Size (above Hip)	NHS Code:
1 x High Receptacle	Small (28-30")	A100S-5105-C
1 x Core Briefs (Variable Size)	Medium (32-34")	A100M-5105-C
1 x Re-useable 500 ml Collection Bag	Large (36-38")	A100L-5105-C
1 x Cleansing Container	XL (40-42")	A100XL-5105-C
1 x Cleansing Solution – 118ml	2X (44-46")	A1002X-5105-C
	3X (48-50")	A1003X-5105-C



Optional Extras: not available on NHS Prescription.



500 mL Bag Sleeve

not available on prescription

One Pack £4.95 + £3 P&P

5 Pack £20* (save £7.95)
(Free delivery)



Secure Wrap

not available on prescription

Pack of 1
£5.95 + £3 P&P

Pack of 4 £20
(save £6.80)
free delivery

iMEDicare Ltd, Unit 11,
Shakespeare Ind. Est.,
Shakespeare Street, Watford, WD24 5RR



EC REP: Eurolink (Europe) Ltd, Avalon House
Marcham Road Abingdon, Oxon OX14 1UD, UK
Manufacturer: Arcus Medical, 2327 Distribution Street
Charlotte, NC 28203 USA

āfex®

is the answer when
Nature Calls..



On NHS FP10 Prescription

Fights Odour, Rash and Embarrassment

Comfortable:
Loose fit, adhesive free receptacle with double-wall design reduces exposure to residual urine.

Discreet:
Briefs conceal system under clothing. Double-sealed collection bag effectively contains odor.

Convenient:
Collection bag connects simply with our twist-n-lock design and easily drains with our one-touch port.

āfex

contact@imedicare.co.uk

+441923 23 7795

iMEDicare

Pelvic Health Naturally...




Adhesive-Free EXternal
Collection System for Men



Post-Prostatectomy Management of Stress Incontinence

www.mypelvichealth.co.uk

Comfortable, Convenient and Cost-effective

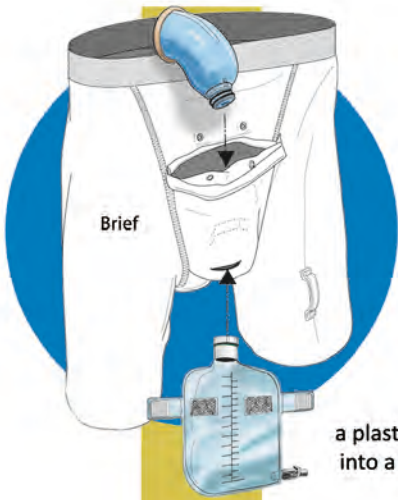

afex[®]

Active / Core

The afex system makes managing post-prostatectomy stress incontinence easier and more dignified than ever. The 96% cotton boxer briefs discreetly and securely hold a plastic (PVC) curved receptacle inside a front pocket that channels urine into a standard 500 mL reusable collection bag (12 ml vinyl - 4 to 6 week lifespan). Uncomfortable tubes or straps are not required - just attached the bag directly to the receptacle through the briefs for hours of use without maintenance or changing !!

The outdoor systems shown here (afex Active & afex Core) can be used for active daytime use such as golfing, fishing, hiking, bowling, walking, gardening, or travelling or any general out-door activity. It has three main components: an ergonomic receptacle, boxer briefs (or jock strap style briefs with afex core) and a collection bag. Receptacles are available in 2 styles (high and low, briefs are based on waist / hip size and collection bags come in 250ml, 500 ml size (Vented for sports, or Unvented as standard) and 1500ml. The annual cost of AFEX to the NHS can be as low as £160 given the durability of components. Free Sample and home assessment available on request.



afex Urine Collection (thigh) Bags (30-40 day lifespan with continuous use)	NHS Code:
250ml Bag	A400S
500 ml Bag	A400-B-08
500 ml Vented Bag	A400-V-08
1500ml Extra Capacity Bag	A400-E-08

Why choose afex?

afex[®] User Success


requires:

- 1) Flaccid Penis to dangle over the testicles when seated
- 2) Reasonable manual dexterity



Pubic Pressure to hold the receptacle against the body wall comes from a tight but comfortable stretch Brief (4% Lycra)

Secure Positioning of the penis centrally within the receptacle by padded secure wrap hugging the shaft of the penis at the open end of the receptacle. Wrap absorbs splash-back and prevents "popping out" with active movement

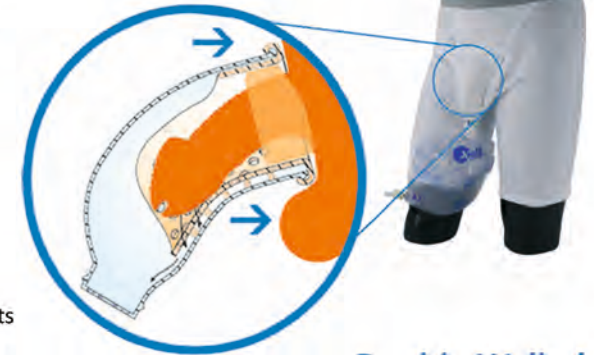

Double Walled

inner receptacle lining drains urine away from the penis, helping to keep skin dry - recovery from sheath excoriation or pad rash.



- ✓ No Adhesives Used
- ✓ Receptacle is Breathable with Natural Relaxed Fit
- ✓ Complete System is Latex Free
- ✓ Simple Three Part System Comes Ready to Use
- ✓ Double-Wall Receptacle Keeps Skin Away From Urine
- ✓ Receptacle with Direct-Connect Bag is Fast Flow Capable
- ✓ Pressure-Free Collection Bag with Activity Tabs

